DENTSPLY

Dentomat[®] compact

Manufacturer:
DeguDent Indústria e Comércio Ltda.
Rua Igarapava, 366
Catanduva - SP
15810-255 - Brazil

EU-Representative:
DENTSPLY DeTrey GmbH
De-Trey-St. 1
D-78467 Konstanz
Tel. +49 (0) 75 31 5 83-0

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ISO 8282:1997

A. Equipment and Parts Supplied

- 1 Dentomat® compact
- 1 Mixing capsule
- 1 Cleaning brush
- 1 Instructions for use
- 1 Certificate of warranty
- 1 Hopper

B. Technical Data

Main power supply: 220 Vac (can be switched to 110 V~)
Frequency: 50/60 Hz
Power consumption: 70 W
Body (h/w/d): 190/145/250 mm
Weight: approx. 4,7kg
Operating temperature: +4 °C - +25 °C
Type: class 4 ISO 8282

C. Description

The Dentomat® compact is a spilling and mixing unit which was developed specifically for the efficient preparation of amalgams. Very fine-grained silver alloys and/or contaminated resp. oxidized mercury can have a negative effect on the proper function. The portions of amalgam can always be dispensed in the correct quantitative ratio, and are homogeneously triturated in preprogrammable times depending on the size.

Design and electrical components meet the requirements for prevention of accidents.

This unit is **not** intended for operation in hazardous areas, i.e., where there is a danger of explosion.

The operating time of the Dentomat® compact must not exceed 1 minute. As a rule, the materials to be processed with this unit require mixing times which lie below the permissible operating time.

We recommend that you unconditionally read the following instructions for use as well as the information relating to maintenance **before the unit is put into initial operation.**

D. Operating Instructions

1. Electrical Connections

The Dentomat® compact is designed for operation of 100 Vac and 220 Vac single-phase, alternating current, frequency 50/60 Hz. The units, when dispatched, are set for 220 V operation. In order to convert a unit to operation at some other voltage, the changeover switch (Fig.1) visible on the bottom of the base plate must be shifted with a screw driver or similar tool. The number that is visible corresponds to the voltage that has been set.

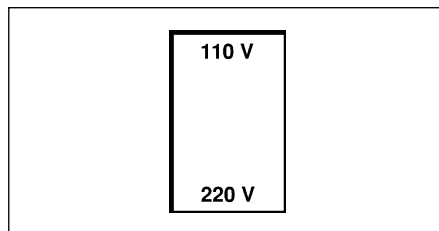
Fig. 1

Switch the instrument on with the "On/Off" switch (left gray button on the frontside), P1 lights up in the digital display (Fig. 2).

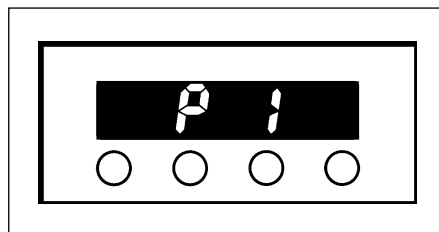
Fig. 2

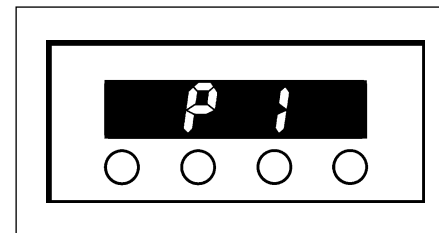
Fig. 2

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2. Filling the Storage Containers

Unscrew the blue screw cover at the top of the unit. Unscrew the seal on the silver alloy container (Alloy) located under the cover. Pour silver alloy powder into the container up to a maximum level even with the lower edge of the white upper section of the dispensing mechanism, and close the container again. The min. filling level is reached when the stirrer is visible

Unscrew the seal on the mercury container, pour mercury into the container. Use a hopper, to avoid the spraying of mercury. Close the container again.

Caution! When filling one supply container, always keep the other container sealed! If the containers are filled with the incorrect material, or if the two charging materials are mixed together, the Dentomat compact will be made unserviceable.

The reservoir may be max. filled to the lowest edge of the cover lid thread. The min. filling level of the mercury reservoir is marked by the top edge of the elongated hole in the reservoir bottom.

Caution statement: Care should always be exercised to prevent mercury spills and leakage. Mercury and mercury dispensers should be kept at less than 25 °C.

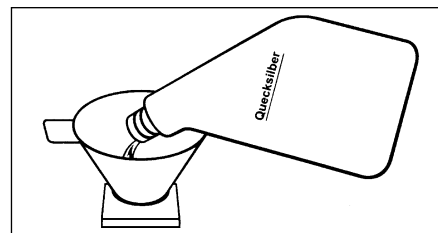
Fig. 3

! Thelf mercury should be spilled during the filling operation, it must be carefully drawn up with an expendable syringe and/or with tin foil (in the case of very small amounts) and stored in a closed vessel under a fixing salt solution! Do not remove it with the saliva ejector!

3. Setting the Mixing Ratio

The mixing ratio of silver alloys powder : mercury is set at the adjustment disc (Fig. 4). For this purpose, the dispensing knob (Fig. 5) must be in the zero position.

Mixing ratios and mixing times are given in the particular working instructions from the producers.

An individual adjustment in the mixing consistency can be made by changing the **amount of powder** ; the mercury portion remains constant and amounts 320 mg each dispensing operation. The adjustment disc (Fig. 4) is turned with a coin.

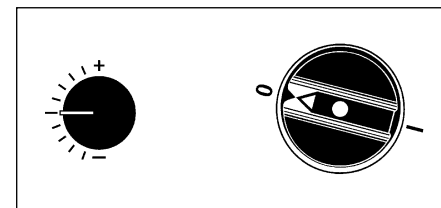
Fig. 4

Fig. 5

Amalgam paste becomes more moist when the mark on the adjustment disc is shifted in the direction of the minus sign.

Amalgam paste becomes drier when the mark on the adjustment disc is shifted in the direction of the plus sign.

! Do not set the mixing consistency too dry. Amalgam paste that is too dry can

not be condensed pore-free even under strong plugging pressure.

The change in the mixing ratio becomes fully effective only after several dispensing operations after the adjustment disc has been reset.


4. Programming the Mixing Times

The Dentomat compact has a storage capacity of 9 different mixing times, P1 - P9.


 The mixing times for the silver alloys you use are given in the particular instructions for use!

4.1 After the Dentomat® compact has been switched on at the "On/Off" switch, the program P1 lights up in the display.


4.2 Press the button "Progr./Time" briefly, and the display shows some time in seconds between "00" and "99".

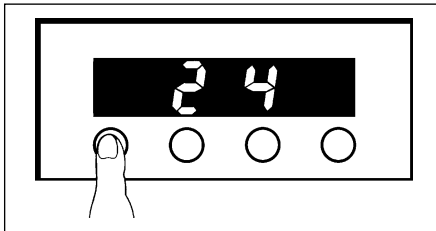
4.3 Enter the 10-second intervals (10, 20, ... 90) with the button "Sec.", for examples.

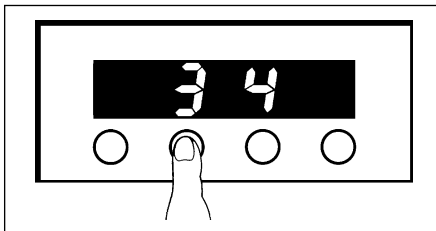4.4 Enter the seconds intervals (1, 2, ... 9) with the button "Progr.-N°", for example 2.

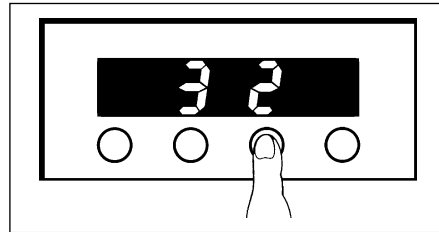4.5 When the button "Progr./Time" is pressed again, "P1" appears in the display, while if the same button is pressed once more, the display shows the time stored in the program 1 (in our example 32 seconds as shown in the figure).

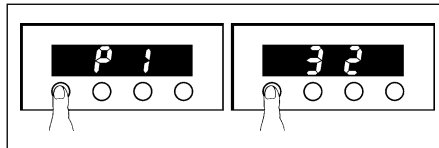4.6 Switch to the next desired program number by pressing the button "Progr.-N°" and proceed as described in sections 4.2 to 4.5.

4.7 After completion of the programming (in general, up to P3 for the usual 3 portion sizes), press the button "Progr./Time" so that P... appears in the display. With the button "Progr.-N°" you can now select the desired program.

Holding the button "Progr.-N°" pressed down for a prolonged length of time causes the program number to switch further, while pressing the button "Sec." causes the program number to switch back.

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
All programmed mixing times remain stored independent of the power supply, i.e., even when the unit is switched off!

When the unit is switched on again by means of the on/off switch, the program P1 always appears, not the last program number that was used last before the unit was switched off.

5. Dispensing

- 5.1 After the Dentomat compact has been switched on, check that the mixing capsule is sealed firmly in place. A mixing capsule that is not screwed firmly in place blocks the unit.


- 5.2 After the selection of a certain program for the spilling and time, operate the button "Start/Stop."

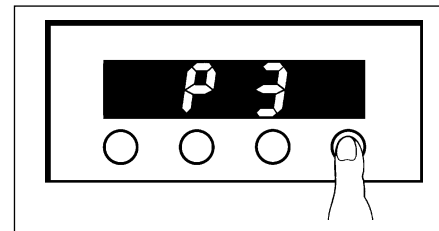- 5.3 The display immediately shows the dispensing time on the program selected, and the mixing capsule vibrates. Now begin immediately with the dispensing.

For this purpose, turn the rotary knob quickly in a clockwise direction as far as the stop marked –, and then turn it back to the starting position 0. As this is done, a distinct snap must be heard at both end points.

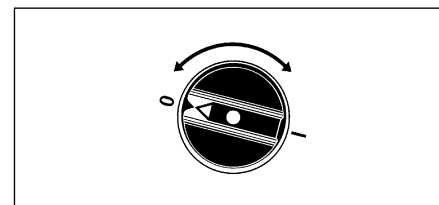Between the two turning movements (initial position 0 => stop – and stop – => initial position 0) there should be a short pause. In order to spill larger portions of amalgam, operate the rotary knob two or three times. Observe the dosing and mixing time.

- 5.4 Removal of the amalgam paste
After the mixing mechanism has been switched off, unscrew the mixing capsule and remove the amalgam paste which is now ready for plugging. Empty the capsule completely. Clean the mixing tube with the cleaning brush, and immediately screw an emptied, clean mixing capsule in place again.

- 5.5 Switch off
Switch the unit off, when it is not to be used for a longtime.
For this purpose push the "On/Off switch once.

If you have inadvertently set the unit in operation by means of the button "Start/Stop", but find that you have selected the incorrect mixing time (the wrong program), interrupt the process by pressing again the "Start/Stop" button. The display shows the mixing time in seconds which is still to be carried out. Now press the button "Progr./Time"; the indication of the previous program then appears in the display. Now, by pressing the button "Progr.-N°" select the correct program.

E. Maintenance and Care


The Dentomat® compact requires occasional maintenance! The mixing tube (connection between the storage container and mixing capsule) has to be cleaned daily with the brush supplied. Never use any probes or similar tools for this purpose because they destroy the high-luster polish of the tube wall.

Mixing capsules which have become worn out, deformed, or contaminated must be replaced by new capsules so that no amalgam residues are mixed in, and in order to maintain the airtight seal of the mixing tube.

The motor and other parts of the Dentomat® compact (including the dosing mechanism) must not be oiled or lubricated! Repairs may be carried out only by the manufacturer or by workshops which are authorized to do this work by the manufacturer.

Before dispatch care should be taken that the machine does not contain any mercury.

⚠ The Dentomat® compact must be protected against extreme fluctuations in temperature (for example, intense sunlight, strong sources of heating).

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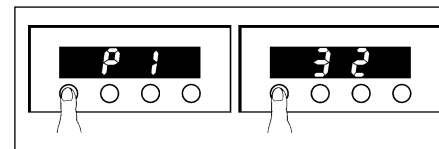This is to certify the conformity of the product

Dentomat® compact / mercury and alloy mixer and dispenser

with the basic requirements of the following EC guidelines:

1. EC guideline for electrical equipment to be used within certain voltage limits (Guideline for Low Voltage) 73/23/EEC with 1st Modification 93/68/EEC
2. EC guideline on electromagnetic compatibility 89/336/EEC with 1st Modification 93/68/EEC

The following harmonized standards were applied:

EN 60204-1:1997; EN 61010-1:2001
EN 61000-1:2001; EN 50082-1:1998

The following international standards and specifications were applied:

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